



Ministère de la
Région Wallonne

ANALYSIS OF THE SEDIMENT FRACTION OF FEED

BY NEAR-INFRARED MICROSCOPY (NIRM)
Vincent Baeten¹, Antoine Michotte Renier¹, Georges Sinaeve¹, Ana Garrido-Varo² and Pierre
Dardenne¹

¹ Agricultural Research Centre (CRA), Quality of agricultural products Department, Chaussée de Namur, 24. B-5030 Gembloux, Belgium. baeten@cra.wallonie.be

² Faculty of Agriculture and Forestry Engineering ETSIAM, University of Cordoba (UCO), Apdo de correos 3019, 14080 Córdoba, Spain.



1. Aim and objectives

With the emergence of the BSE crisis we assist to the development of a wide range of analytical methods to detect meat and bone meal in feed ingredients and compound feeds (Gizzi *et al.*, 2003). Indeed, contaminated feeds are commonly accepted as the main transmission of the mad cow disease in the European bovine herds. Briefly, there are the methods based on optical microscopy, molecular biology (e.g. PCR and ELISA) and spectroscopic analysis of the sample. Classical microscopy analysis of the feed is the reference method used throughout Europe to control the correct implementation of the ban of meat and bone meal. Firstly, it was applied for the ruminants and then extended to all farmed animals. With this method, the detection of meat and bone meal is done by microscopic observations, after adequate sample preparations including the use of staining reagents, of specific features of particles coming from different sieved and/or sediment fractions. The test is mainly based on the detection of bones by visual observation of the particles resulting from the sedimentation of the sample.

The main spectroscopic methods proposed to control the illegal addition or to detect the contamination of feeding stuffs are based on the near-infrared spectroscopic technique (Garrido and Fernandez, 2000; Baeten and Dardenne, 2001). With these methods, the infrared spectra from the raw sample (near-infrared spectroscopy) or from the particles making up the sample (near-infrared microscopy) are used to detect samples including meat and bone meal material. In previous studies (Piroux and Dardenne, 1999; Baeten *et al.*, 2001), it has been demonstrated the high potential of near-infrared microscopy method to detect meat and bone meal at level as low as 0.5%. To decrease the limit of detection and to increase the speed of the near-infrared microscopy analysis, the addition of a sedimentation step before the near-infrared microscopic analysis has been investigated.

2. Results and discussion

In order to calibrate the near-infrared microscope for the detection of bones in the sedimented fraction, a series of feed ingredients from animal, vegetal and mineral origin as well as MBM free compound feed were sedimented according to the procedure described in the UE microscopic guideline and were analysed by near-infrared microscopy (NIRM).

Figure 1 presents several microscopic images and NIR spectra of particles coming from the sedimented fraction of a compound feed. The animal particles spectra shows bands with maximum near 1946, 2044, 2176 and 2288 nm.

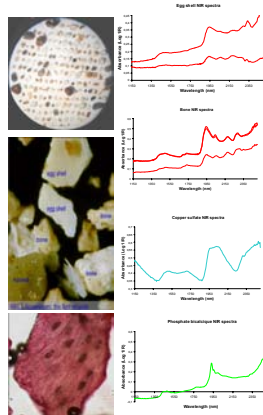


Figure 1. Microscopic images and infrared spectra of different particles of the sedimented fraction. (Part of the microscopic images are issued from the Decision Support System ARIES develop in the framework of the STRATFEED project by RIKILT (Dr. Leo van Raamsdonck) and aiming to assist the microscopist in the detection and identification of animal protein in feeding stuff)

2. Results and discussion (cont.)

To test the method developed at CRA, a spiked sample at 0,1 % with MBM was analysed ten times. **Table 1** presents the results. The number of analysed particles, the number of particles detected as having an animal origin, the percentage of bone particles in the sediment and the weight of bones in the sediment are included. All the replicates were detected as positive regarding the presence of bones in the sample. The standard deviation of the percentage of bone particles in the sediment (SD = 0.16) underlines the quality of the method proposed taking into account the difficulty to get a perfect homogeneous spiked sample at this low level.

Table 1. Analytical results of the 10 NIRM analysis of a sample spiked at 0.1 % of meat and bone meal.

Code	% sediment (a)	Nb analysed particles	Nb animal particles (b)	% bones in the sediment (c)	Weight bones in the sample (d) [mg]	Conclusion
a	2,73	333	2	0,63	1,73	Positive
b	2,57	329	3	0,96	2,47	Positive
c	2,59	340	2	0,62	1,61	Positive
d	2,96	330	2	0,64	1,89	Positive
e	2,81	305	2	0,69	1,94	Positive
f	2,66	292	2	0,72	1,92	Positive
g	2,55	342	1	0,31	0,79	Positive
h	2,49	329	2	0,64	1,59	Positive
i	2,72	279	2	0,75	2,05	Positive
j	2,37	378	2	0,56	1,32	Positive

Legend: (a) % sediment = (Weight sediment/Weight sample) * 100
(b) Nb animal particles = Number of particles clearly identified (by means of their infrared spectrum) as having an animal origin
(c) % bones in the sediment = (Nb animal particles/Nb analysed particles) * 100; % of particles of the animal food ingredients detected as animal origin
(d) Bones in the sample = Weight sediment * % bones in the sediment (c)

In the framework of the STRATFEED project (<http://stratfeed.cra.wallonie.be>), a total of 24 MBM free compound feeds and 48 spiked samples at a level ranging from 0.5 to 8 % were analysed using the NIRM method. The samples were prepared by UCO and analysed in blind by CRA. All the free MBM compound feeds were detected as negative samples and all the adulterated samples were detected as positive regarding the presence of meat and bone meal. **Figure 2** displays the analytical results of a series of the blind samples spiked with the same MBM.

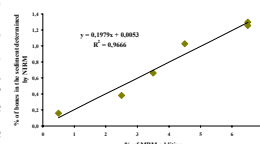


Figure 2. NIRM results of the analysis of a series of blind samples.

4. Conclusions

The analysis of the sedimented fraction by near-infrared microscopy is an attractive and powerful method to detect low level of meat and bone meal in feeding stuffs. The main features of this method are:

- the analysis of particles from the sample. This is similar to the classical microscopy (UE reference method).
- the limit of detection (LOD < 0,1 %).
- the use of infrared detectors and discriminant equations instead of the eyes and the expertise of the analyst as within the classical microscopic method.

Reference

- Gizzi, G., van Raamsdonck, L.W.D., Baeten, V., Murray, I., Berben, G., Brambilla, G., von Holst, C. [2003]. An overview of tests for animal tissues in animal feeds used in the public health response against BSE. *Sci. Tech. Rep. In press*.
- Baeten, V. & Dardenne, P. (2001). The contribution of near infrared spectroscopy to the fight against the mad cow epidemic. *NIRS news*, 12 (6), pp. 12-13.
- Baeten, V., Michotte Renier, A., Sinaeve, G. & Dardenne, P. (2001). Analysis of feeding stuffs by near-infrared microscopy (NIRM): detection and quantification of meat and bone meal (MBM). *Proceeding of the 6th International Symposium on Food Authenticity and Safety FADIS conference*, Nantes, France, November 28-30 2001, pp. 1-11.
- Garrido-Varo A. (2000). - La spectroscopie infrarouge: une technologie d'appui pour un service integral en alimentation animale. *In La spectroscopie infrarouge et ses applications analytique* (D. Bertrand & E. Dufour, eds), éditions TEC & DOC, Paris, Chapter 20, 473-496.
- Piroux F. & Dardenne P. (1999). - Feed authentication by near-infrared microscopy. *In Proceeding of the 9th International Conference on Near Infrared Spectroscopy* (A.M.C. and Giangiacomo R. eds), Verona, Italy, June 1999.

Acknowledgement

Part of this project is financed by the European Commission, UE STRATFEED project **G6RD-2000-CT-00414**. Strategies and Methods to detect and quantify mammalian tissues in feedingstuffs. <http://stratfeed.cra.wallonie.be>.